

Applications of TCR repertoire analysis for biomarker discovery and beyond

Ankita Das¹, Vasumathi Kode¹, Kayla Lee¹, Priyanka Shah², Xiaoshan "Shirley" Shi¹, Nitin Mandloi, Ravi Gupta², Amit Chaudhuri^{1,2} and Papia Chakraborty¹

1) MedGenome Inc, Foster City, CA, 2) MedGenome Labs, Bangalore, India; corresponding author : ankita.d@medgenome.com

Abstract

T cell immunity provides significant therapeutic benefit to cancer patients treated with checkpoint inhibitors, however a very small fraction of patients typically respond to checkpoint inhibitors, and a smaller fraction of them have any long-term benefit. This can be attributed to the lack of prognostic and predictive biomarkers. The infiltration levels of CD8 T cell in tumors is often used as a characteristic biomarker and can be correlated with response, but recent studies have shown that examining the functional state of the T cells and other immune cell types in the tumor, and the immunogenic neoantigen burden and the TCR repertoire clonality, might give a more appropriate representation of what might be going on in the tumors and hence can be used as predictive biomarkers for response. At MedGenome we have built a suite of tools that can be utilized to: a) Study the tumor microenvironment using the OncoPeptTUME analysis pipeline, b) Predict and validate the immunogenic neoantigens using OncoPeptVAC & OncoPeptSCRN, and c) A suite of workflows to analyze the TCR repertoire from a wide range of sample types using bulk and single cell approaches. Here we present the data highlighting the applications of TCR repertoire as a biomarker for immunotherapy and also present our capabilities of providing these solutions as a service.

Overview of the tumor extrinsic and intrinsic factors that determine immunotherapy success

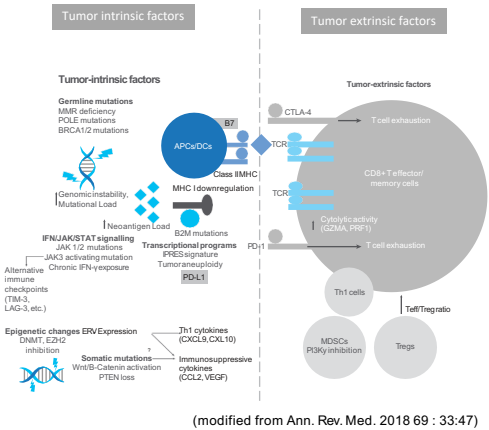


Figure 1: Shows an overview of the complex interplay between the tumor cell and the immune cell types that can be in the tumor, and the cell extrinsic and intrinsic factors in a patient's tumor that can be important in determining success to immunotherapy or not. On the left there are some of the tumor intrinsic molecular features: for instance if a tumor carries certain germline mutations that result in genomic instability, leading to high tumor mutational load and some of those mutations leading to expression of neoantigens in the tumor, if those neoantigens are then presented on the surface of the cell by HLA presentation, and presented via an APC, then cells in the microenvironment such as a CD8 T cell can be attracted to that target and induce tumor killing. However, they need to have the appropriate TCR repertoire to recognize the HLA bound peptide antigens to induce downstream signaling and mount an anti-tumor response. There are other immune cell types in the tumor such as MDSCs and T-regs that can be suppressive for inducing tumor killing mechanisms. Therefore, by studying the immune infiltration in the tumor, and identifying potential immunogenic neoantigens, and looking at the TCR repertoire can all be informative biomarkers for immunotherapy.

MedGenome's tools to enable to biomarker discovery for Immunotherapy

Table 1: Summary of the biomarker discovery services at MedGenome

Biomarker/therapeutic	MedGenome's solution	Deliverables
Tumor microenvironment analysis	OncoPeptTUME	Gene expression analysis, immune score, cytokines & pathways
TCR repertoire analysis	OncoPeptTCR	V(D)J gene clusters, CDR3 clonotypes, TCR α/β pairing
Neoantigen prediction and validation	OncoPeptVAC & OncoPeptSCRN	Predicted & validated neoantigen vaccine candidates

Overview of the workflow for OncoPeptVAC pipeline and TCR repertoire analysis

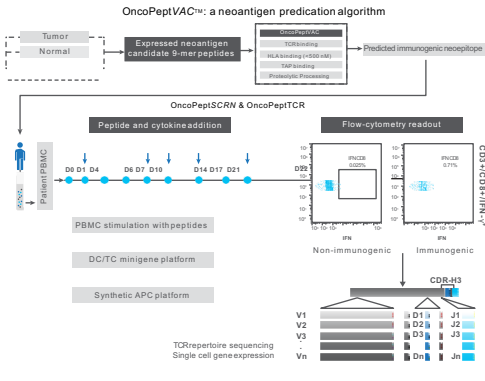


Figure 2: Workflow of neoantigen prediction and screening: The OncoPeptVAC pipeline takes whole exome sequencing and RNA-sequencing data from a performed using the genomic DNA from blood and tumor and RNA from the patient's tumor. All mutant peptides are passed through the OncoPeptVAC pipeline. The potential immunogenic peptides are predicted based on a series of prioritization steps including peptide processing, peptide TAP binding, HLA binding and TCR binding. Finalized peptides are then subjected to OncoPeptSCRN wherein peptides are tested for CD8 T cell activation assay in presence of antigen loaded APCs. Functional TCR cell response is assessed by IFN γ production by flow cytometry, TCR repertoire sequencing and single cell gene expression analysis. The TCR repertoire and single cell gene expression analysis are performed to reveal the clonality and the functional T-cell states of the amplified T-cells.

TCR repertoire sequencing combined with single cell gene expression analysis reveals the clonal amplification and functional T-cell states

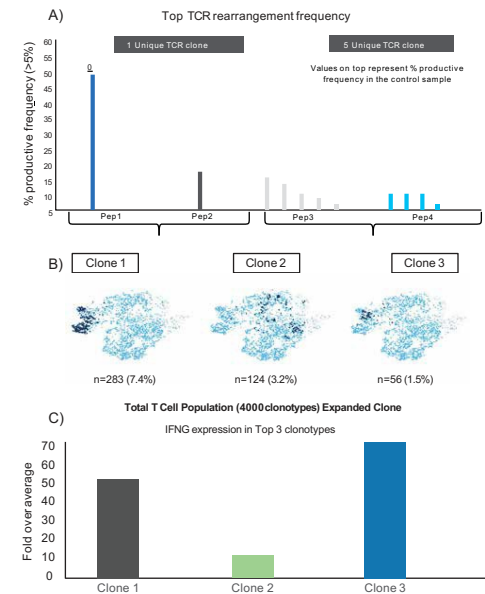


Figure 3: TCR repertoire and single cell gene expression analysis to identify antigen specific TCR repertoire clonal expansion: A) TCR repertoire analysis of PBMCs screened using the neoantigen vaccine candidates shows that compared to their wildtype counterparts, mutant forms of the peptide candidates show TCR α/β clonotype frequency changes between 20%-50%. We also observed that while some vaccine candidates mount a monoclonal response, others induce a polyclonal response. The TCR repertoire libraries were generated using the SMARTer TCR α/β Profiling Kit (Takara Bio USA), and sequencing was performed using the Illumina MiSeq platform, and repertoire analysis was performed using the MixCR pipeline. B&C) Single cell TCR repertoire analysis of the PBMCs (with positive responses to the peptide, assayed by IFN γ production) identified over 4000 TCR clonotypes, with the t-SNE plots showing the top 3 clonotypes and the frequencies are listed below. Single cell gene expression of the amplified T cells reveal that there is varying levels of IFN γ in T-cells that are clonally expanded, and that this information can provide insights on the functional state of the expanded T-cell clones. The single cell TCR and gene expression experiment was performed using the 10x Genomics Immune profiling solution & sequencing was performed using the Illumina NovaSeq 6000, and 10x Cell Ranger analysis was performed downstream.

γ/δ TCR repertoire Profiling as a biomarker & potential therapeutic

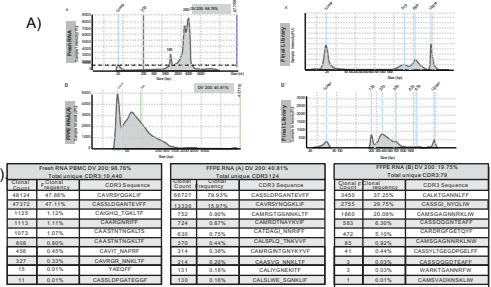


Figure 4: TCR repertoire analysis from FFPE tissue: A) Tapestation profiles of the RNA from PBMC and an FFPE sample and the corresponding TCR α/β libraries generated from each of the samples. The expected TCR libraries from good quality RNA is 500-700bp, and with FFPE a skewing of the peaks between 200-500bp is observed. B) Table shows the top clonotypes identified from the three different FFPE samples. We show that for all the samples regardless of DV200, we identify over a 100 clonotypes. In addition to that, we observed that multiple samples from patients treated with the same vaccine show a shared clonotype among the top clonotypes observed. The FFPE TCR assay was performed using 10ng of total RNA, and a workflow developed by MedGenome using the SMARTer TCR α/β Profiling Kit (Takara Bio USA), and sequencing was performed using the Illumina MiSeq and analysis was performed using the MixCR pipeline (MLabs).

γ/δ TCR repertoire Profiling as a biomarker & potential therapeutic

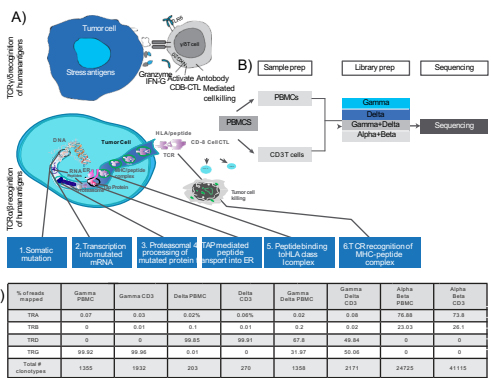


Figure 5: TCR repertoire analysis from γ/δ repertoire: A) Cartoon representation of the γ/δ TCR compared to the α/β TCR, where in contrast to MHC-presentation of antigens, γ/δ T cells are able to recognize stress antigens without MHC-presentation, making them attractive for cancer immunotherapy. They can have a CTL phenotype via engagement with CD8 T cells. B) Shows workflow for enriched specificity of γ/δ TCR assay using 10ng of total RNA from PBMCs and CD3 enriched T cells. C) Table shows proof of concept data showing that in TCR libraries where γ/δ primers were added we see specific enrichment of γ/δ chains versus alpha beta and similarly in TCR α/β libraries we see specific amplification of TCR α/β chains. This assay is performed using 10ng of total RNA from PBMCs and CD3 sorted cells, and γ/δ primers designed and optimized by MedGenome using the SMARTer Takara TCR α/β Profiling Kit, and sequencing was performed using the Illumina MiSeq, and analysis was done using MixCR.

Summary of the TCR sequencing services at MedGenome

Table 4: Summary of the for TCR repertoire sequencing at MedGenome

Technology	Sample & Input Type	Reagents	Sequencing Method	Analysis Method	Information Obtained	Amount of Input Material	Approximate Sequencing Cost	Read Length (bp)	Number of Clonotypes
Single Cell TCR	Individual Cells (PBMCs)	SMARTer TCR α/β Profiling Kit	Illumina MiSeq	MixCR	CD8α/β pairing and clonotype V(D)J sequences	~100,000 cells	\$500	150bp	~100,000 clonotypes
Bulk SMARTer TCR Profiling	Isolated Cells (PBMCs)	SMARTer TCR α/β Profiling Kit	Illumina MiSeq	MixCR	CD8α/β pairing and clonotype V(D)J sequences	~100,000 cells	\$500	150bp	~100,000 clonotypes
FFPE Bulk TCR Profiling	FFPE RNA (FFPE)	SMARTer TCR α/β Profiling Kit	Illumina MiSeq	MixCR	CD8α/β pairing and clonotype V(D)J sequences	~100,000 cells	\$500	150bp	~100,000 clonotypes
Single Cell BCR	Individual Cells (PBMCs)	SMARTer BCR α/β Profiling Kit	Illumina MiSeq	MixCR	CD8α/β pairing and clonotype V(D)J sequences	~100,000 cells	\$500	150bp	~100,000 clonotypes
Bulk SMARTer BCR Profiling	Isolated Cells (PBMCs)	SMARTer BCR α/β Profiling Kit	Illumina MiSeq	MixCR	CD8α/β pairing and clonotype V(D)J sequences	~100,000 cells	\$500	150bp	~100,000 clonotypes

Table 5: Summary of the sample types for TCR

Sample Type	Approximate Sequencing Cost	Read Length (bp)	Number of Clonotypes
Single Cell TCR	\$500	150bp	~100,000 clonotypes
Bulk SMARTer TCR Profiling	\$500	150bp	~100,000 clonotypes
FFPE Bulk TCR Profiling	\$500	150bp	~100,000 clonotypes
Single Cell BCR	\$500	150bp	~100,000 clonotypes
Bulk SMARTer BCR Profiling	\$500	150bp	~100,000 clonotypes